


Changing vest liner on Halo

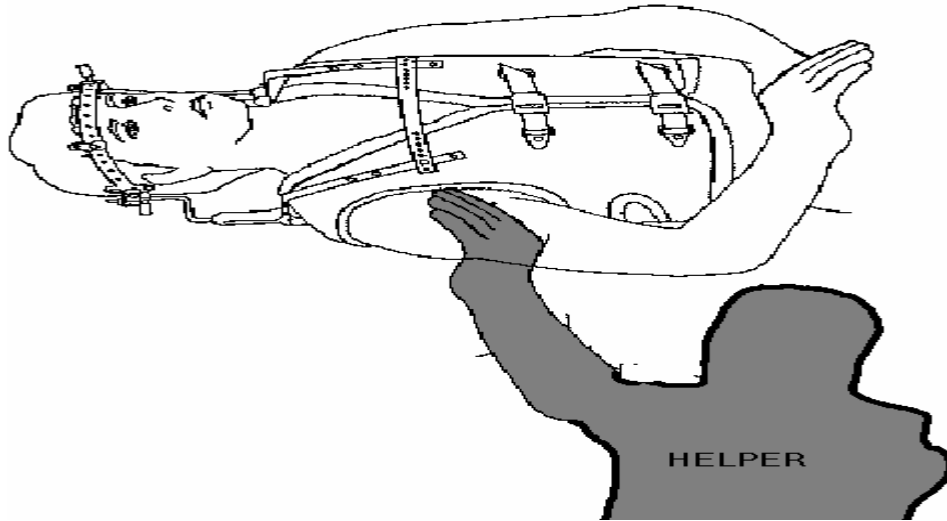

With person lying on their back:

Slide your hands in between the plastic part of the vest and the sheepskin.

You can work your hand around and pull the sheepskin down from the velcro attached to the inside of the plastic

As you work your hand around you should be able to separate the layers.

To add the new layer of sheepskin. It is best to slide a pillow case or towel, or plastic bag in between the plastic and the body.

After the pillow case is in place slide the sheepskin between the pillow case and the body. (the pillow case keeps the velcro from “gripping” before you have the sheepskin in place)

After the sheepskin is in place pull the pillowcase out and press the sheepskin into the velcro

To do the back sheepskin piece it is possible to have the person lie on their side.

Questions:

Call MOPC @ 800 598 6672